

The following guide provides a brief overview regarding assembly and disposal of 5-quart sharp containers provided by Hazardous Materials Management.




Assembly and Use

Press the lid into the base making sure all six tabs are fully seated into the base and have snapped fully into the base. The six tabs should be seen completely seated through the underside of the red base (pictured left). Please do not fill above the fill line (pictured left) as this may hinder the closure of the rotating lid.

Disposal

- Submit the sharp container at utsa.SciShield.com, using the sharp container profile name.
- Verify “rotating cylinder lid” tabs (2) are snapped into the lid.
- Place a piece of tape over the lid to further latch the sharp container (pictured below)
- Write the request ID# from SciShield on the sharp container (pictured below).

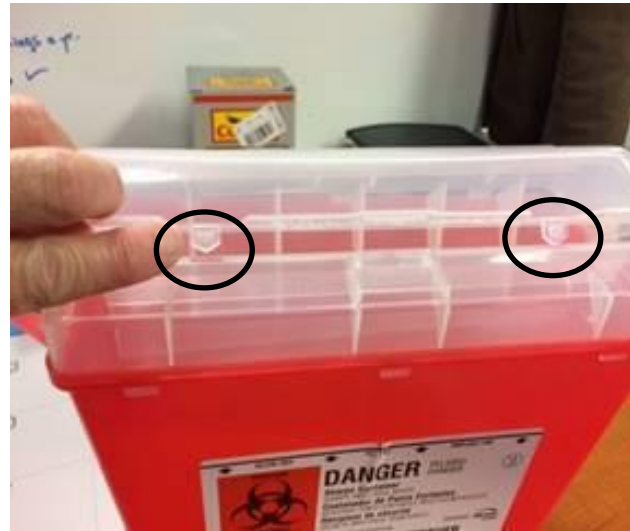

Notes

- Sharp containers are given out on a one-for-one basis
- No liquid-containing ampules, test tubes, syringes, etc should be placed in the sharp containers. All liquid should be collected in another container and submitted as chemical waste.
- Failure to follow disposal procedures may result in rejection of waste submission.